
THUASNE

SLEEQ[®] MAX


- (EN) Premium Post-Op/Post-Injury TLSO – INSTRUCTIONS FOR USE
- (DE) Premium TLSO nach Operationen/Verletzungen – Gebrauchsanweisung
- (FR) Corset thoraco-lombo-sacré (TLSO) post-opératoire/après blessure haut de gamme – Mode d'emploi
- (IT) Istruzioni per l'uso – Ortesi premium di stabilizzazione toraco lombo sacrale post-opertoria/post lesione
- (ES) Premium Post-Op/Post-Injury TLSO – Instrucciones para su uso
- (PT) TLSO premium pós-operação/pós-lesão – Instruções de utilização
- (NL) Premium TLSO voor post-operatieve toepassingen/behandeling van letsel – Gebruiksaanwijzing
- (DA) Kvalitets-TLSO efter operation/skade – Brugsanvisning
- (SV) Premium post-op/post-skada TLSO – Bruksanvisningar
- (FI) Leikkauksen/vammanjälkeinen selkärankatuki – Käyttöohjeet
- (EL) Premium μετεγχειρητική/μετατραυματική θωρακο-οσφυο-ιερά όρθωση (TLSO) – Οδηγίες χρήσης
- (CS) Pooperační a poúrazová ortéza TLSO – Návod k použití
- (PL) Pooperacyjna/pourazowa orteza krzyżowo-lędźwiowo-piersiowa (TLSO) Premium – Instrukcja użytkowania
- (ET) Kvaliteetne operatsioonijärgne / vigastusjärgne TLSO – kasutusjuhend
- (LV) Augstākās kvalitātes pēc operācijas/pēc traumas TLSO – Lietošanas norādījumi
- (LT) Aukščiausios rūšies pooperacinis / potrauminis TLSO (torakolumbosakralinis ortezas). Naudojimo instrukcijos
- (SK) Prémiová pooperačná/poúrazová torako-lumbo-sakrálna ortéza – Návod na použitie
- (SL) Visokokakovostni pripomoček TLSO za postoperativno uporabo in uporabo po poškodbah – navodila za uporabo
- (HU) Prémium műtét/sérülés utáni TLSO – Használati útmutató
- (RO) Orteză toraco-lombo-sacrală (corset TLSO) premium postoperatorie/post-traumatism – Instrucțiuni de utilizare
- (BG) Висококачествени тораколумбални сакрални ортези (TLSO) за след операции/след наранявания – Инструкции за употреба
- (HR) Vrhunska TLSO ortoza za postoperativno razdoblje ili razdoblje poslije ozljede – Upute za upotrebu
- (ZH) 优质的术后/受伤后胸腰骶护腰 (TLSO) – 使用说明

2.0 SIZING

- 2.1 PREPARE TO FIT:** Push button (1)
Remove anterior panel (A).

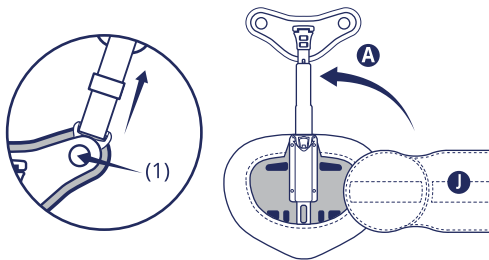
- 2.2 ADJUST CIRCUMFERENCE:** Measure waist and see Sizing Guide (1.0).

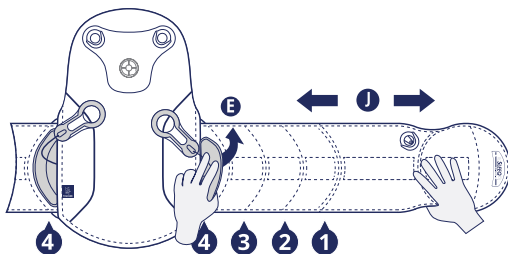
3.0 ATTACH WING SHELLS & LATERAL PANELS


- 3.1 ALIGN** wing shell (C) sew line Arcs to sew line Arcs on wing. **Note:** For size (4) use posterior Arc near Velcro tab (A).

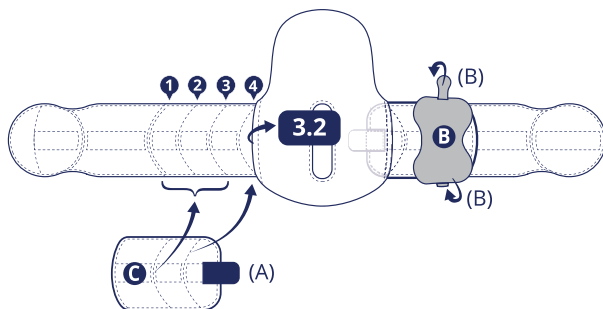- 3.2 TUCK** wing shell (C) inside back panel.

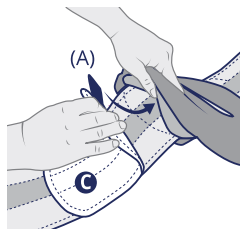- 3.3 SECURE** tab (A) to wing inside back panel. - Repeat on opposite side.

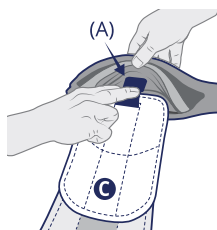
- 3.4 ALIGN** Lateral Panel (B). Attach by folding over both tabs (B).


4.0 ADJUST AND FIT

- 4.1 POSITION** brace at coccyx level.

Note: Have patient hold wings.

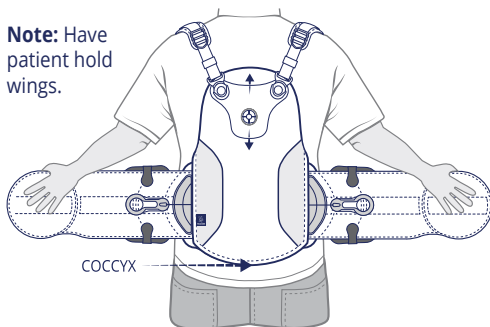
- 4.2 POSITION** Anterior Panel (A). Wrap wing (1) and attach to anterior panel. Wrap wing (2) and secure.

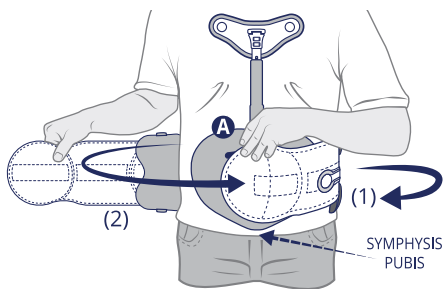
4.0 ADJUST AND FIT (Continued)

4.3 TIGHTEN BRACE: Move pull-tabs **D** anteriorly and re-attach.

4.4 STERNAL BAR: Open sternal bar clip **L** to adjust height, then close clip to lock.
Sternal plate **K** is 1 to 2" (25-50cm) below Sternal Notch when standing and 1/2" (23cm) below when seated.
Sternal bar **M** can be bent to achieve proper fit with bending irons or over the edge of a table.
Note: Only bend Sternal Bar above toothed area.

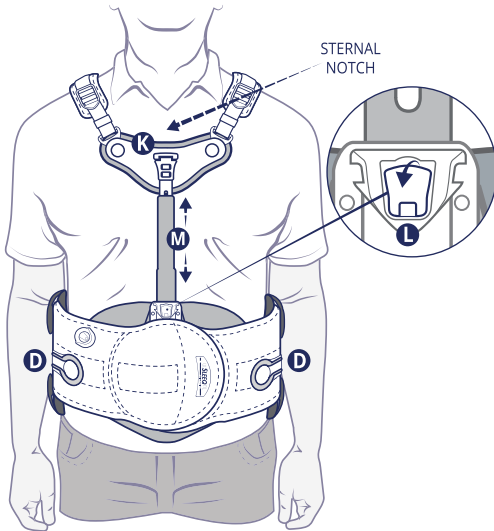4.3 REMOVING OR REPLACING STERNAL BAR (with available Short Bar option)

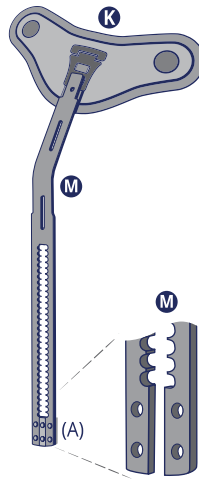- Remove anterior panel from wing.
- Release sternal bar clip **L** and lower bar.
- Remove retainer clip **(A)** by pushing through bar while bending bar away from panel.
- Pull bar **M** up to remove from bar slot.
- Replace bar, or (optional Short Bar) into bar slot and replace retainer clip **(A)**.


Note:
Clip **(A)** must be replaced on the panel side of bar.

5.0 Modifications

The Back Panel, Anterior Panel **A** (shown), Lateral Panels and Wing Shells can be trimmed with scissors **(1)** and heat modified **(2)** to adjust lordosis, add scapula reliefs or create flares. The steel Sternal Bar **M** can be bent to match anatomical contours and the wings can be attached to the Velcro® Connectors at up to 5 degree angles up or down to accommodate various waist to hip ratios.

CAUTION: Do not trim top center tabs **(3)** on anterior panel.

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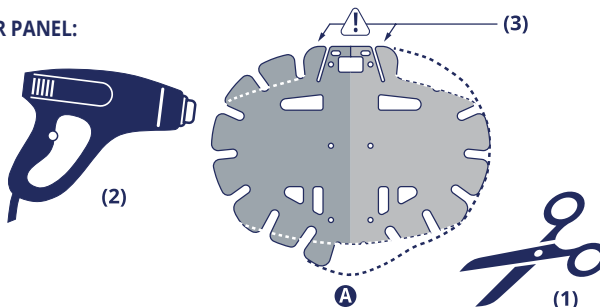
5.0 Modifications (Continued)


- 5.2

BACK PANEL REMOVAL FOR MODIFICATION:
(1) Open cover at the base of the back panel. (2) Lift cover. (3) Remove the LSO assembly. (4) Remove plastic panel.


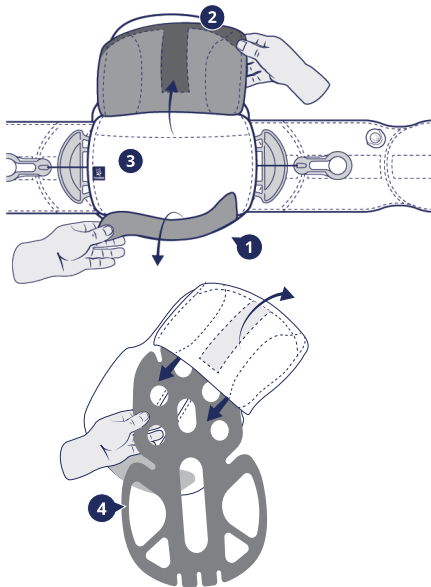
- 5.3

ADJUST FOR OVER THE SHOULDER OR UNDER AXILLA STRAPPING OPTION
Detach Strap Anchor **N** and reposition up for over the shoulder (**A**) or down for under axilla strap placement (**B**).

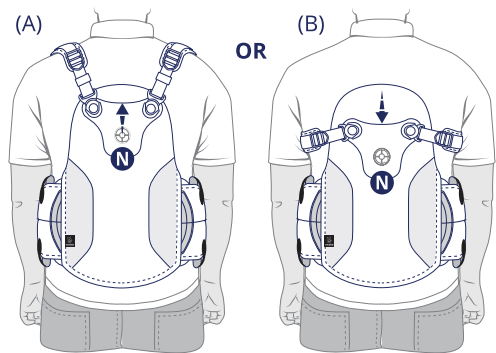6.0 Reset

When removing the brace, always replace Pull Tabs **D** to ready position (6.1). Lay brace on a flat surface then pull both Wings **I** outward to reset the tightening mechanism (6.2).

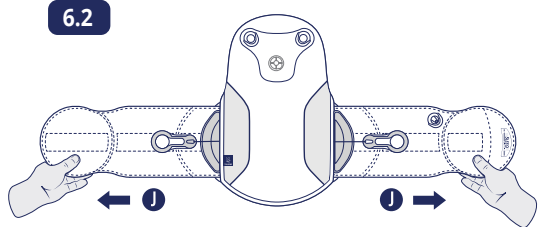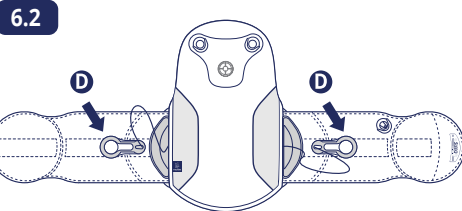SIDE EFFECTS:

This device can cause skin reactions (redness, itching, burns, blisters...) or wounds of various degrees of severity.

WARNINGS:

- Carefully read use/care instructions and warnings prior to use.
- If increased pain, swelling, sensation changes, or any adverse reactions are experienced while using this product, immediately consult your medical professional.

PRECAUTIONS:

Follow the advice of the professional who prescribed or supplied the product. Consult your professional if discomfort occurs. Store at room temperature, preferably in its original packaging. For hygiene and performance reasons, do not re-use the product for another patient. Any serious incident that occurs in connection with the product must be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient are/is domiciled.

CARE:

Do not use any cleaning agents, conditioners or aggressive products (containing chlorine).

Squeeze out the water thoroughly.

Dry away from any heat sources.

**MATERIAL COMPOSITION:**

Nylon, Polyester, Polypropylene, synthetic rubber (latex-free), neoprene, stainless steel, brass, aluminum.

DISPOSAL:

Dispose of in accordance with local regulations.

KEEP THIS INSTRUCTION LEAFLET.**Limited Warranty:**

Thuasne USA will repair or replace all or part of the unit and its accessories for material or workmanship defects for a period of six months from the date of sale.

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