



THUASNE

BOLD
LIGAMENT BRACING

 **TOWNSEND**
THUASNE USA

BOLD

Patients who have experienced ligament injuries and/or reconstructive knee surgery.

FEATURES

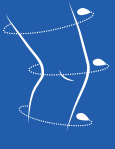
- Lightweight brushed aluminum frame with black polymer cuffs
- Townsend Motion TM5+ Hinges
- Exclusive Synergistic Suspension Strap
- 14 Inch length
- Extension stop kit
- Anti-migration strap padding

OPTIONS

- Flexion stop kit
- Knee guard

ORDERING

- OTS: SM to XL

Size	7" above	7" below	
SM	15.5" - 18.5"	12.25" - 13.75"	
MD	18.5" - 21"	13.25" - 15"	
LG	21" - 23.5"	14.25" - 15.75"	
XL	23.5" - 25"	15" - 17"	



The BOLD Difference...

SELF-FORMABLE POLYMER THIGH AND TIBIA CUFFS

The semi-malleable upper and lower cuffs conform intimately to the leg – achieving enhanced support and patient comfort.

EMBOSSED AIRCRAFT ALUMINUM SUPERSTRUCTURE

Embossing the rigid shells substantially increases the BOLD's structural integrity for contact sports and active patients.

TOWNSEND HINGE TECHNOLOGY

The BOLD integrates Townsend's roll-and-glide hinge motion. Tracking with the arthrokinematic motion of the knee provides continuous brace-to-leg contact throughout range of motion, enhancing control and functional mobility.

SUSPENSION TECHNOLOGY

In concert with the hinge mechanics, Townsend's patented Synergistic Suspension Strap prevents brace migration and enhances patient satisfaction. Strategic attachment and angulation of the strap creates circumferential compression and utilizes the calf muscle as a natural shelf.

MOISTURE WICKING PAD LINERS

Like the tread on a car tire, the pad design features channels that allow moisture to exit the liner – enhancing suspension and patient comfort.

