



**THUASNE**

## Ligaflex Boxer and Ligaflex Finger Braces







# Ligaflex® Boxer

## Wrist-Hand-Finger Orthosis

PDAC Approved L3807/L3809

### Indications

- Non-displaced 4th and 5th metacarpal fracture
- Postoperative 4th and 5th metacarpal support
- Following strict cast immobilization

### Comfort

- Soft fabric with foam liner to cover the injured area
- 3D inner breathable fabric for daily comfortable wear
- Finger fabric separator to protect the immobilized finger

### Immobilization

- Two shapable rigid stays for a strict immobilization from the wrist up to the finger to aid in healing
- Pre-shaped palmar and dorsal stays (dorsal = 90° and palmar = 20°) for optimum immobilization

### User-Friendly

- Full opening of the brace for easier application and positioning
- Convenient buckle loop system to tighten the brace
- Machine washable



# Ligaflex® Finger

## Wrist-Finger Orthosis

Recommended PDAC L3923/L3924

### Indications

- Long finger injuries including sprains, dislocations and fractures (functional treatment or post surgical support)
- After-effects of injuries (sprain, dislocation and fracture)

### Comfort

- Four straps, two fixed and two removable, for optimized support and fit to improve wearing comfort
- Wrist strap with anatomical shape to improve wearing comfort
- Curved top of the stay to cover and protect the finger

### Immobilization

- Unique stay with four straps to stabilize the brace
- Rigidity for strict immobilization from the wrist up to the finger
- Adaptability to fit the patient's wrist and hand shape and follow the finger angulation
- Two stay positions, palmar and dorsal, for optimum immobilization

### User-Friendly

- Full opening of the brace for easier application and positioning
- Machine washable








## Ligaflex® Boxer

Wrist-Hand-Finger Orthosis


Size	Wrist Circumference
1	5 – 5 3/4"
2	6 – 6 3/4"
3	7 – 7 3/4"
4	8 – 9"

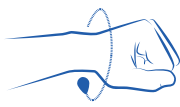

## Ligaflex® Finger

Wrist-Finger Orthosis

Size	Tip of Injured Finger
1	< 5"
2	5 – 6 1/4"
3	6 1/4 – 7 1/2"
4	7 1/2 – 8 3/4"

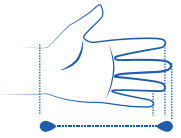

**Townsend /Thuasne USA**

ThuasneUSA.com • 4615 Shepard Street, Bakersfield, CA 93313 • 800.432.3466