



THUASNE

---

# ROM LX

PDAC L1832/1833





Full Pads Available  
(Q3 2022)


# ROM-LX

## PDAC APPROVED L1832/1833

**Recommended for:** Patients who require post-operative fixation and/or range of motion control for post-injury stabilization.

### Feature & Benefits

- ROM hinges with tool-free adjustable drop locks that can be engaged at all increments
- Simple flexion and extension pull pin system that engages in 10° increments
- Aluminum struts with ease of use telescoping ends that allow simple height adjustment ranging from 17 ¼" – 26 ½" in 1/16" increments.
- Malleable paddles with individual cool pads that can be contoured to the leg to achieve an intimate fit, increase stabilization, and improve brace suspension
- Universal size, no measurement required

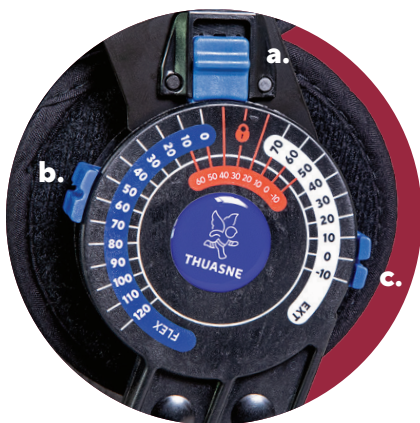

#### Drop Lock (a)

0° to 120° in  
10° increments

#### Flexion Stop Setting (b)

0° to 120° in  
10° increments

#### Extension Stop Setting (c)

-10° to 60° in  
10° increments



Button Control  
Length Adjustment



**Thuasne USA**

4615 Shepard Street, Bakersfield, CA 93313 • 800.432.3466

[ThuasneUSA.com](http://ThuasneUSA.com)