



THUASNE

SLEEQ[®] MAX

Post-op/post-injury TLSO
with universal sizing




Recommended for: Moderate to severe sagittal and coronal transverse plane motion restriction / increased intracavitary pressure to reduce spine and disc loading / mechanical or discogenic thoracic lumbar and extremity pain / post-operative spinal stabilization / fusion / degenerative disc disease / post-injury to facilitate healing / spondylolysis / spinal stenosis / spondylolisthesis / spondylosis, laminectomy / discectomy / facet syndrome / spinal muscle weakness / arthritis / low back sprains, strains and prophylactic back support.

SLEEQ[®] MAX PDAC Approved L0464

Rigid tri-planar control

Achieve maximum control with the SLEEQ[®] Max's versatile strapping system, sternal bar, full circumference shell design and reinforced panels.

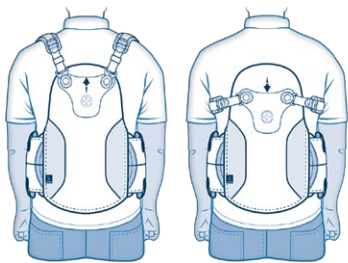

Over the shoulder or under axilla strapping options

Modular step-down design

From TLSO to LSO, support your patient throughout the recovery process. Removable straps & back panel, and a detachable sternal bar & pectoral plate means the SLEEQ Max can transition from a L0456 to a L0627 in minutes.



Adjustable compression

Smooth action pull tabs adjust the corset style lacing system and allows patients to easily control the amount of compression and support.

Comfortable and easy to wear

Increase patient compliance with the TLSO that's simple to don and wear. Featuring an ergonomic strap system, self molding plastics, and breathable material, the Max makes wearing a TLSO easy!



Simple operation of snap and release buckles

Universal sizing with InvisAdjust[™]

Fits 26"-54" waist circumference and up to 62" waist with available EXT extension panel.



Steel sternal bar assembly

- Pressure relief sternal pad
- Steel upright with snap-lock height adjustment
- Short sternal bar option for smaller patients
- Self molding steel reinforced anterior panel




Modular rigid back panel

- Repositionable harness with posterior and anterior strap clips
- Rib reinforced anatomically contoured panel with 15° lordosis
- Under axilla strapping option
- Steps down to LSO




Under axilla straps maximize flexion control

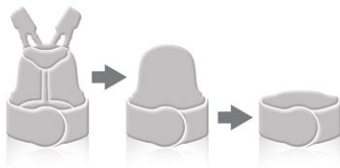


Full-circumference shell design

- Reinforced wings and lateral panels create continuous support
- Dressing relief cutout

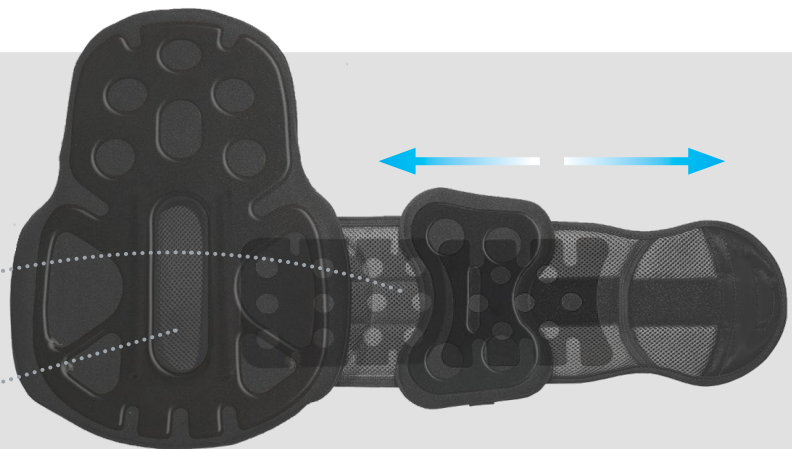

Customizable components

Smart, modular design for an array of modifications. Trim, bend, or heat mold plastic rigid interior panels to adjust lordosis, increase pressure relief, or create flares. Steel sternal bar can be bent to match anatomical contours. Wings can be angled up or down to support wide hips or a pendulous abdomen.

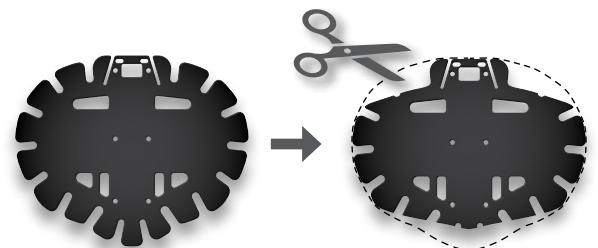

Anterior panel can be trimmed with scissors

Partner Program

Your logo on every SLEEQ® Max brace. Ask your Thuasne USA representative about the free program that puts a spotlight on your brand.



(((SLEEQ)))

SPINAL BRACES

A comprehensive range of Universal LSOs and TLSOs that provide versatility, comfort, and ease-of-use. All SLEEQ® models are customizable by practitioners and offer built-in compression adjustment for maximum patient comfort.

More options. Better care and support.



MODEL	SLEEQ AP	SLEEQ AP+	SLEEQ APL	SLEEQ FLEX	SLEEQ MAX
PART #	901100	901200	901300	901400	901600
PDAC	L0627/L0642	L0631/L0648	L0637/L0650	L0456/L0457	L0464

4615 Shepard Street, Bakersfield, CA 93313
800.432.3466 • Fax: 800.798.2722

ThuasneUSA.com