



THUASNE®

MOBIDERM® STANDARD



*Move with Comfort,
Dream in Style*

The **MOBIDERM**[®] Standard

The Mobiderm[®] Standard garments are expertly crafted to target and mobilize edematous or hardened / fibrotic subcutaneous tissue - ideal for the maintenance phase of lymphedema treatment once limb volume has stabilized or plateaued. Designed with patient compliance in mind, these garments offer a soft, comfortable touch in a slim, modernized design. Simple to slip on styles to encourage autonomy. Let your patients rest as the garment gets to work!



Phase of Treatment: Maintenance



Time of Day: Evening, Nighttime, and Low Activity



Manufacturing Type: Off-The-Shelf



Color: Beige




Warranty: 6 Month Manufacturer's Warranty



The **MOBIDERM**[®] Technology

Mobiderm[®] foam blocks provide a comfortable, massaging effect to support clinical efficacy and compliance. This technology is made up of big or small foam blocks encased between two non-woven bandages. Due to the difference in pressure created between the support area of the blocks and their surrounding area, it promotes the flow of lymphatic fluid and optimizes the drainage efficacy. ⁽²⁾

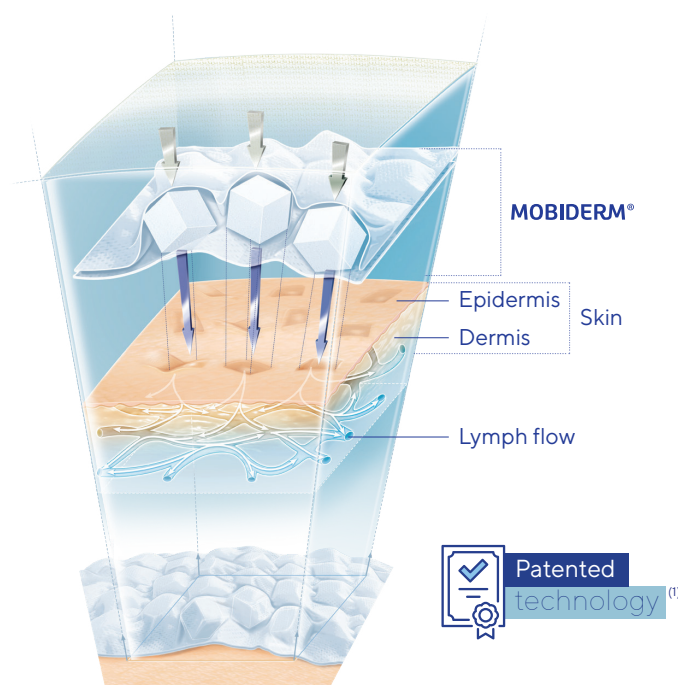

What to Love About Mobiderm® Standard?

Proven Efficacy of Mobiderm® Technology

- ✓ Reinforce the effects of the daytime treatment ⁽²⁾⁽³⁾
- ✓ Designed to maintain the volume reduction of the edema ⁽³⁾ and lymphedema
- ✓ Promotes the flow of lymphatic fluid and optimizes the drainage efficacy ⁽²⁾
- ✓ Mobilizes edematous or subcutaneous tissue (e.g. fibrosis)



Results after just one night

PATIENT TESTIMONIAL | Tina Jackson

"Finally after 15+ years living with lymphedema from breast cancer, I found the perfect night garment. The way it hugs my arm and compresses it so well has made a significant difference in reducing the tightness and thickness of my swelling. It is light weight compared to my other night garments and I am excited to use this as my new travel garment as well."

THERAPIST TESTIMONIAL | Barbara Tobin, OTR/L, CLT-LANA (Daughter of Tina)

"As a lymphedema therapist, it's essential to me that my patients feel both comfortable and confident in their night garments. Through my own experience with my mom, I've seen firsthand how challenging living with lymphedema can be. But I've also witnessed the powerful impact of positive reinforcement and the right garments—ones that fit well, feel comfortable, and truly make a difference. The clear improvement in my mom's arm, visible in the garment's indentations that it made, speaks volumes. I'm truly grateful for Thuasne."

Mitten



Ambidextrous

Turns inside out for right or left usage in one design

Mobiderm® Small Block

Small foam blocks are better adapted to smaller anatomies, such as the dorsum of the hand, top of the wrist and circumferentially around the fingers

Discreet

Slim and flexible design achieved by having Mobiderm® Small Blocks on the dorsum side of the hand and circumferentially around the fingers (glove only)

Extended Wrist


Provides great coverage to be worn alone or overlap with the Mobiderm® Standard Armsleeve

Glove



Separated Fingers

Only nighttime garment with individual, separated fingers supporting excellent mobility

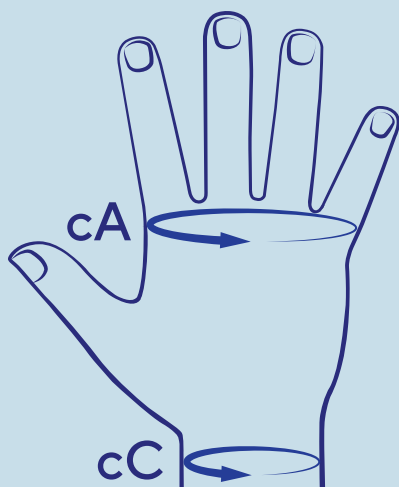

cA: Hand circumference
(head of metacarpal bones)

cC: Wrist circumference



Mitten

Size	cC (in cm)	cA (in cm)	Reference
			right/left
1	14-16	15-17	373101 20100299
2	16-18	17-19	373101 20200299
3	17-19	19-21	373101 20300299
4	19-21	21-23	373101 20400299
5	20-22	23-25	373101 20500299
6	21-23	25-27	373101 20600299



Glove

Size	cC (in cm)	cA (in cm)	Reference	
			right	left
1	14-16	15-17	373201 20100201	373201 20100202
2	16-18	17-19	373201 20200201	373201 20200202
3	17-19	19-21	373201 20300201	373201 20300202
4	19-21	21-23	373201 20400201	373201 20400202
5	20-22	23-25	373201 20500201	373201 20500202
6	21-23	25-27	373201 20600201	373201 20600202

Armsleeve



Mobiderm® Big Block

Featuring big foam blocks for addressing fibrosis and improving lymph drainage

Great Flexibility

Garment is slim and flexible allowing greater movement

Extended Sizes


Includes regular and plus sizes for the upper arm

Slim Design

Discreet treatment in a slim, low profile garment



Armsleeve

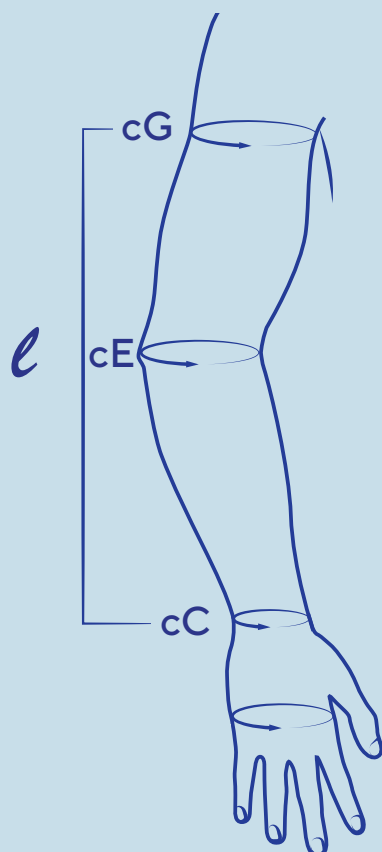

Size	cC (in cm)	cE (in cm)	cG (in cm)	Reference (right/left)	
				Normal ℓ = 40 - 45 cm	Long ℓ = 45 - 50 cm
1	14-16	22-25	25-29	373002 20100232	373002 20100233
1+		24-27	28-32	373002 40100232	373002 40100233
2	16-18	25-28	29-33	373002 20200232	373002 20200233
2+		27-30	31-35	373002 40200232	373002 40200233
3	18-20	28-31	32-36	373002 20300232	373002 20300233
3+		30-33	35-39	373002 40300232	373002 40300233
4	20-22	31-34	35-39	373002 20400232	373002 20400233
4+		33-36	38-42	373002 40400232	373002 40400233
5	22-24	34-37	38-42	373002 20500232	373002 20500233
5+		36-39	41-45	373002 40500232	373002 40500233

cC: Wrist circumference

cE: Elbow circumference

cG: Arm circumference

ℓ: Length from wrist to axilla



Care Instructions

It is recommended to wash your Mobiderm® garment daily (or after each use).

- Use mild detergents. Avoid use of bleach, fabric softener or other harsh chemicals.
- Hand wash:
 - Soak in cold soapy water
 - Rub with a clean sponge
- Machine wash:
 - Gentle / Wool Cycle on cold setting
 - Place the garment in a wash bag
- Dry:
 - Squeeze out excess water
 - Air dry, away from direct heat (e.g. radiator, sun, etc.)
 - Do not tumble dry
- Do not iron
- Do not dry clean



Materials

Outer Layer:

- Polyamide
- Elastane

Inner Layer

- Mobiderm® Big Block:
 - Adhesive non-woven, polyurethane foam
- Mobiderm® Small Block:
 - Adhesive non-woven, polyethylene foam



Indications

- Lymphedema
- Edema

Wings for your health!



Scan Here
to learn more about the
Mobiderm® Standard and
access ordering info.

Thuasne, LLC.

120 Osage Ave.
Kansas City, KS 66105
United States

Phone 1.800.821.3094

Fax 1.800.462.4707

Email cs@therafirm.com

ThuasneUSA.com